AccuCARE™

REFERENCE GUIDE

Equipment:

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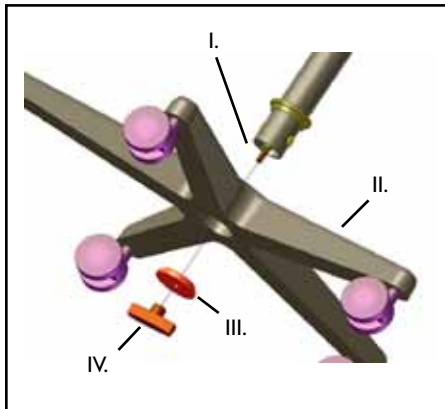
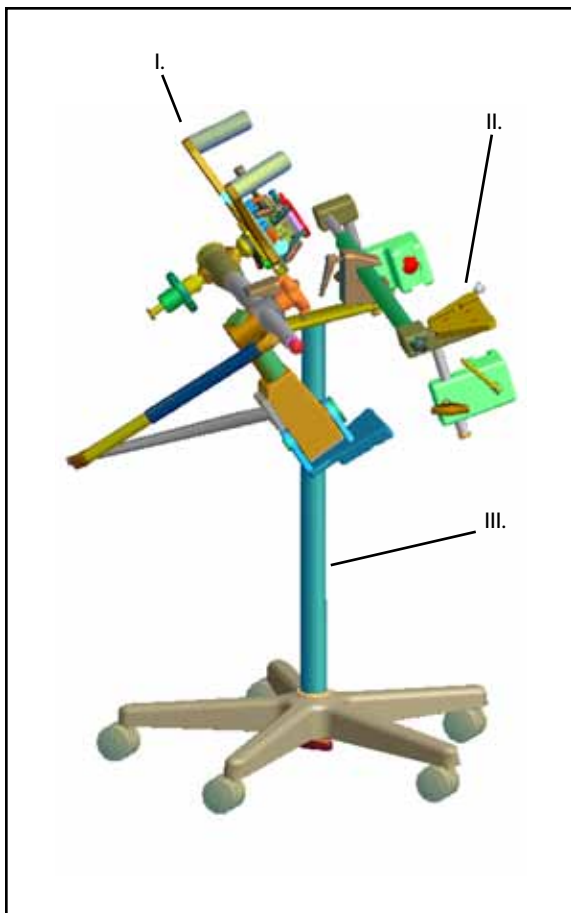
ASSEMBLE TRANSPORTATION STAND**For use with Micro-Touch (Dual-sided table mount)**

1. Insert threaded stud end of tube through center of rolling base.
2. Slide locating washer over threaded stud and nest into center of hole in rolling base.
3. Turn locking knob onto threaded stud and secure firmly.

- I. Threaded stud
- II. Rolling base
- III. Locating washer
- IV. Locking knob

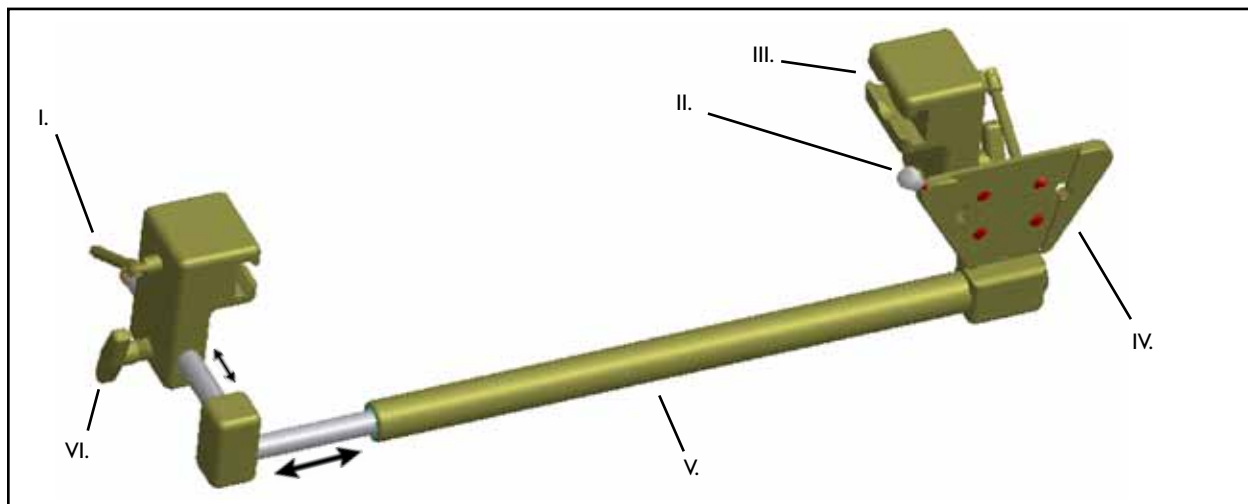
4. The transportation stand provides a convenient and mobile platform to allow ease of cleaning and handling within the surgical environment. Stand shown with stabilizer and table mount clamped in place.

- I. Stabilizer
- II. Table mount
- III. Transportation stand

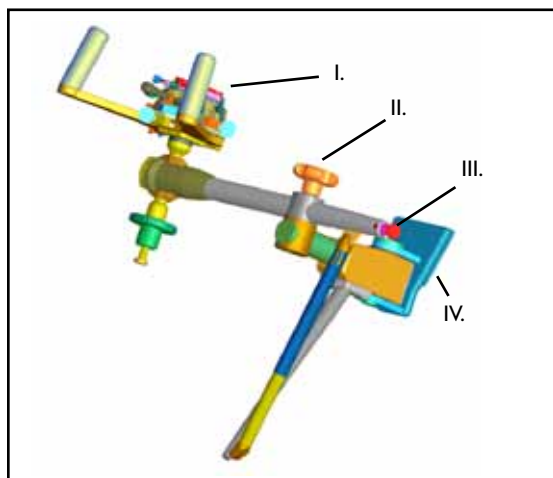
SETTING UP MICRO-TOUCH (DUAL-SIDED TABLE MOUNT)

- I. Loosen all knobs and levers on table mount before removing from transportation stand. Grasp rail clamps and slide entire apparatus into position on end of table. Lock rail clamps onto table by tightening rail clamp levers on both sides at similar distance from the end of rails. Position crossbar at desired distance from end of table to avoid any interference with table mechanisms and tighten both crossbar positioning knobs to secure.

**WARNING**

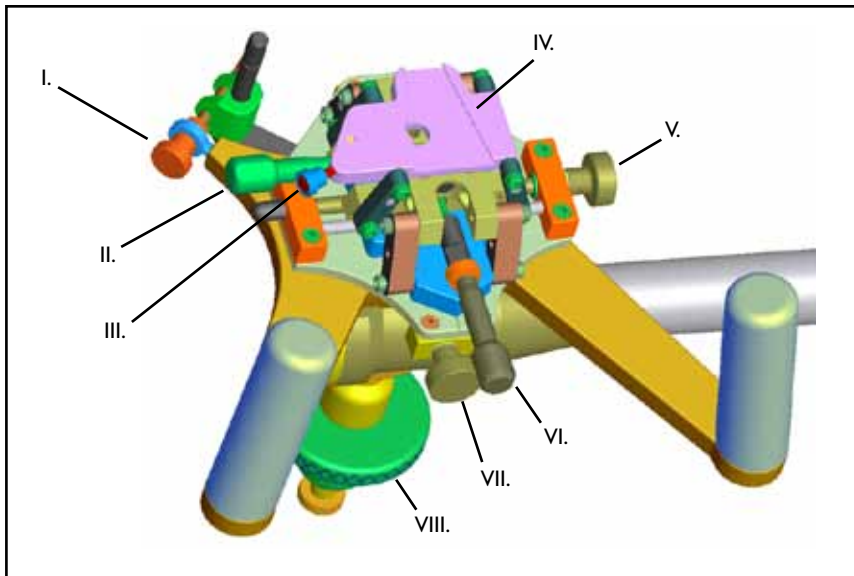
- This system consists of crossbar with attached rail clamps and dovetailed mounting plate for stabilizer. The rail clamps are designed to fit the rectangular side-rails of standard sized operating and cystoscopy tables.

- I. Rail clamp levers
 - II. Stabilizer quick connect knob
 - III. Rail clamp
 - IV. Stabilizer mounting plate
 - V. Crossbar
 - VI. Crossbar positioning knobs
2. Remove stabilizer from transportation stand. Fit dovetailed mounting plate to mating stabilizer mounting plate on table mount. Secure with stabilizer quick connect knob. Release locking control knob and manually position the fine adjustment mechanism (FAM) in a convenient position to facilitate positioning of fine adjustment controls and attachment of stepper. Tighten locking control knob to secure.

- I. Fine Adjustment Mechanism (FAM)
- II. Locking control knob
- III. Weight adjustment knob
- IV. Dovetail mounting plate

3. Check to see all FAM controls are set near midpoints of their range of travel.

- I. Knob YR - Rotates probe to parallel with prostate centerline
- II. Knob ZR - Levels prostate horizontally on transverse image
- III. Stepper Quick Connect Knob - Secures stepper to FAM
- IV. Stabilizer interface plate
- V. Knob X - Centers probe along midline axis of the prostate and centers the transverse image on grid
- VI. Knob XR - Tilts probe parallel to the posterior surface of the prostate in the sagittal plane
- VII. Knob Z - Used to "zero" the stepping unit at the bladder neck to set the base plane for implantation
- VIII. Knob Y - Adjusts image height on grid to define posterior row for seed placement

NOTE | • Knob Y controls vertical movement and has a useful range of about 6cm (3 up and 3 down). Going down it will bottom out as lower plates come into contact and it is shipped in this position. To setup for use, turn Knob Y until lower plates are separated by 2-3cm.

FUNCTIONAL CHECKS

PRIOR TO USE

- Attach the table mount at a suitable clasp point on the table rails to ensure it remains securely in position.
- Ensure stabilizer moves freely when locking control knob is loosened and secures firmly when locking control knob is tightened.
- Move each of the knobs on the FAM through its full range and then center each knob in position.

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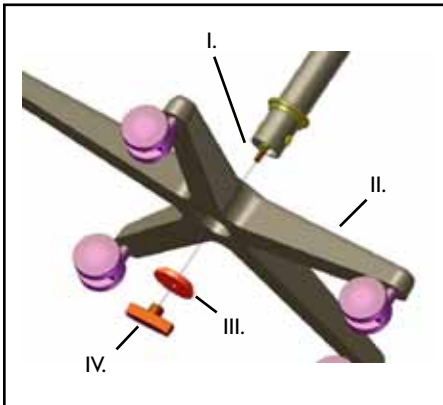
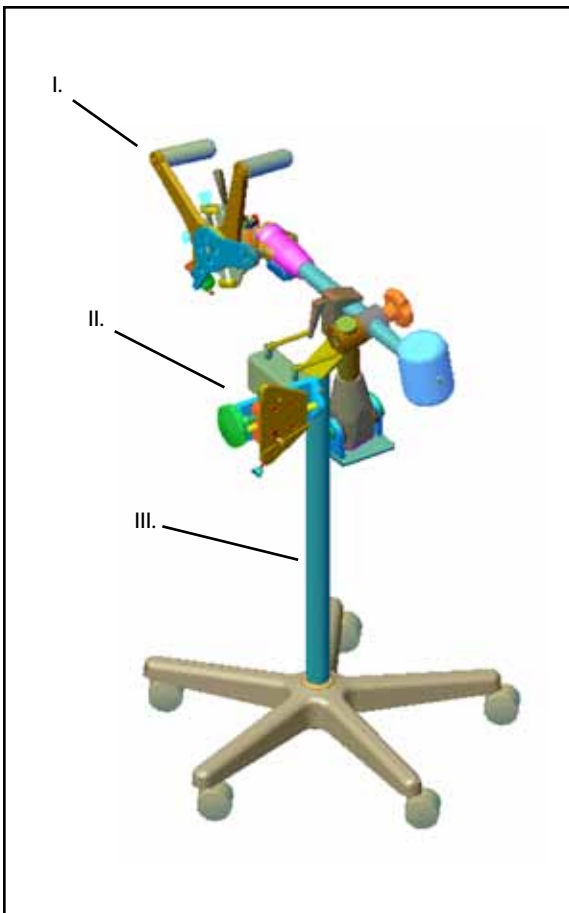
ASSEMBLE TRANSPORTATION STAND**For use with Micro-Touch® LP (Single-sided table mount)**

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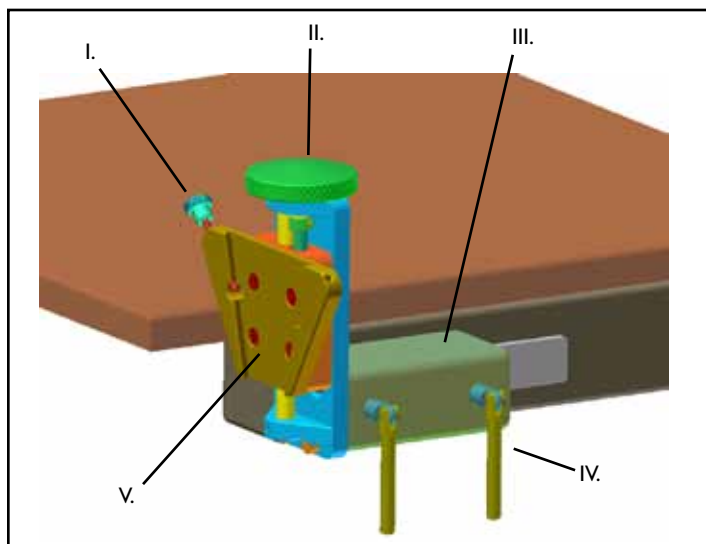
- I. Stabilizer
- II. Table mount
- III. Transportation stand

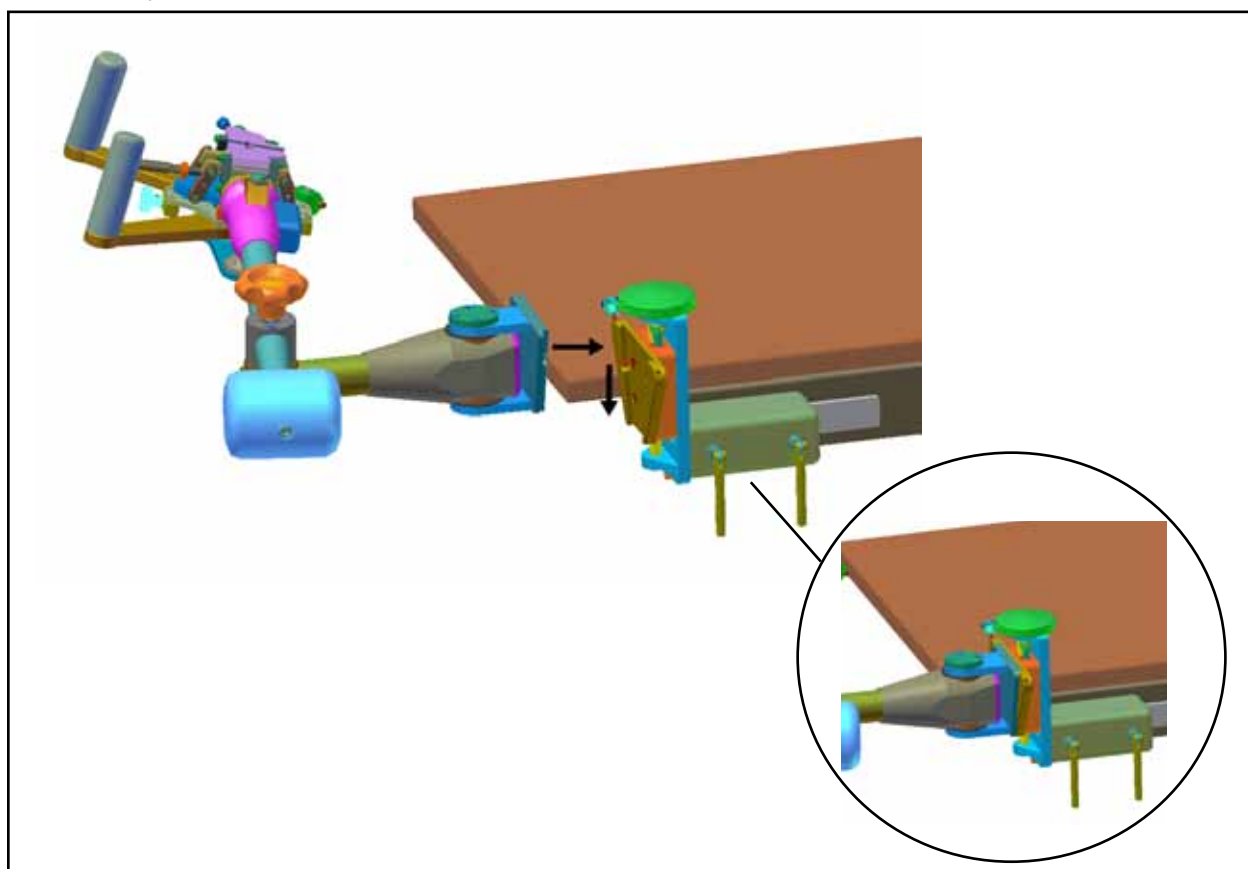
SETTING UP MICRO-TOUCH® LP (SINGLE-SIDED TABLE MOUNT)

- I. Loosen rail clamp levers on table mount before removing from transportation stand. Slide table mount onto table rail and tighten rail clamp levers to secure. Center the Y adjustment.

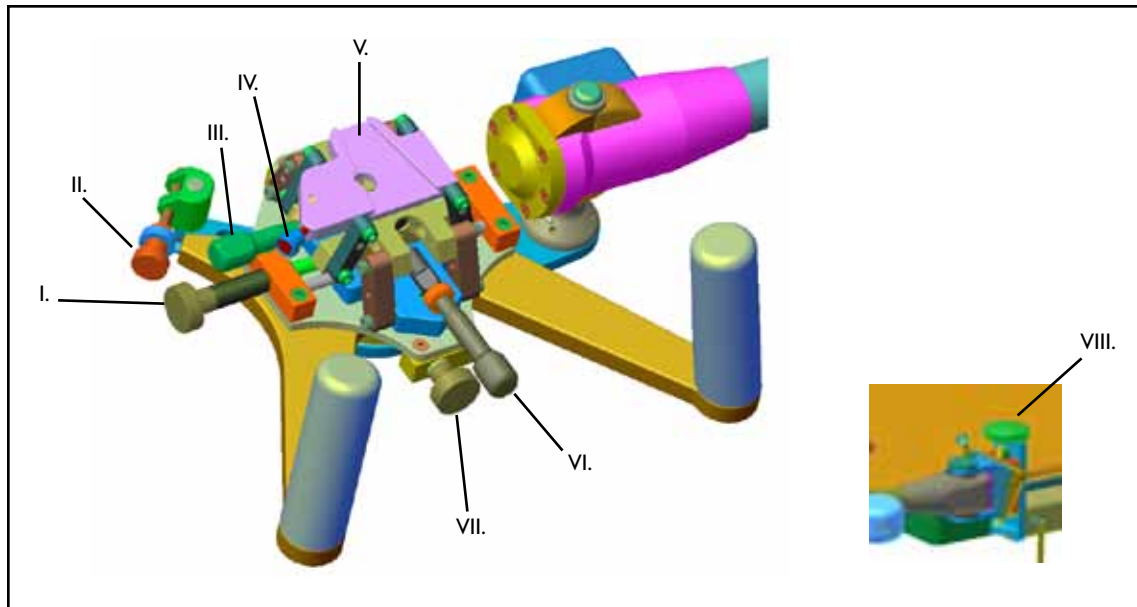
- I. Stabilizer quick connect knob
- II. Knob Y
- III. Rail clamp
- IV. Rail clamp lever
- V. Stabilizer mounting plate

2. Remove stabilizer from transportation stand. Fit dovetailed mounting plate to mating stabilizer mounting plate on table mount. Secure with stabilizer quick connect knob.

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- I. Knob X - Centers probe along midline axis of the prostate and centers the transverse image on grid
- II. Knob YR - Rotates probe to parallel with prostate centerline
- III. Knob ZR - Levels prostate horizontally on transverse image
- IV. Stepper Quick Connect Knob - Secures stepper to FAM
- V. Stabilizer interface plate
- VI. Knob XR - Tilts probe parallel to the posterior surface of the prostate in the sagittal plane
- VII. Knob Z - Used to "zero" the stepping unit at the bladder neck to set the base plane for implantation
- VIII. Knob Y - Adjusts image height on grid to define posterior row for seed placement (Located on the outboard end of the rail clamp.)

NOTE • Due to low clearance of sliding tabletop in relation to fixed components in some tables, a foam mattress pad may be required to allow complete freedom of table movement when the single-sided table mount is in operating position with ultrasound transducer in rectum.

FUNCTIONAL CHECKS

- PRIOR TO USE**
- Attach the table mount at a suitable clamping point on the table rails to ensure it remains securely in position.
 - Ensure stabilizer moves freely when locking knob is loosened and secures firmly when locking knob is tightened.
 - Move each of the knobs on the FAM through its full range and then center each knob in position.
 - Perform table clearance check. Before indexing the table, make sure stabilizer and table mount do not interfere.

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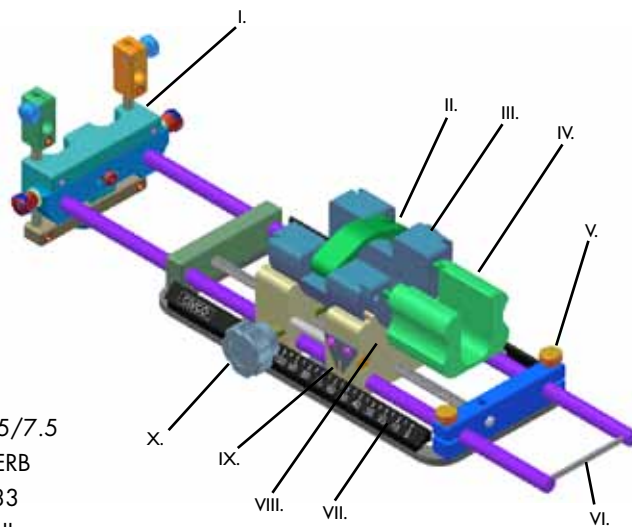
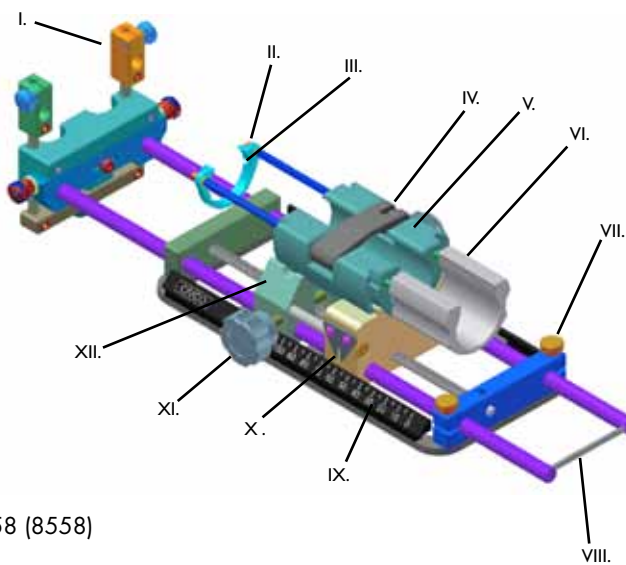
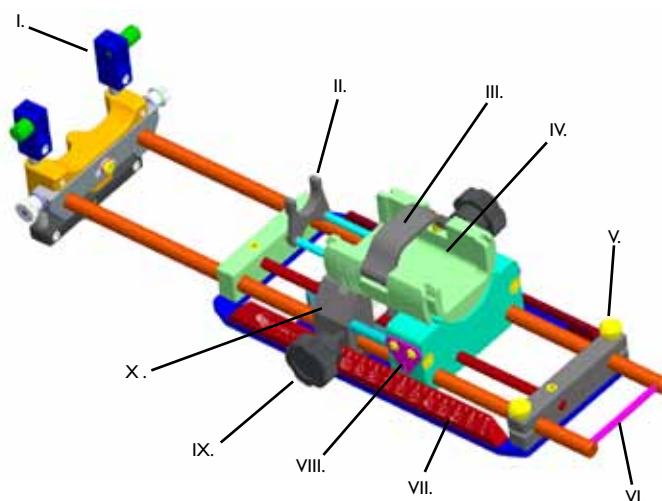
- I. Grid platform
- II. Cradle latch
- III. Cradle
- IV. Cradle handle (*not included on all models*)
- V. Grid rail locking knobs
- VI. Grid rail handle
- VII. Carriage scale
- VIII. Carriage
- IX. Carriage scale marker
- X. Carriage knob

B-K Medical 8658 (8558)

- I. Grid platform
- II. Probe positioning locator
- III. Probe positioning notch
- IV. Cradle latch
- V. Cradle
- VI. Cradle handle
- VII. Grid rail locking knobs
- VIII. Grid rail handle
- IX. Carriage scale
- X. Carriage scale marker
- XI. Carriage knob
- XII. Carriage

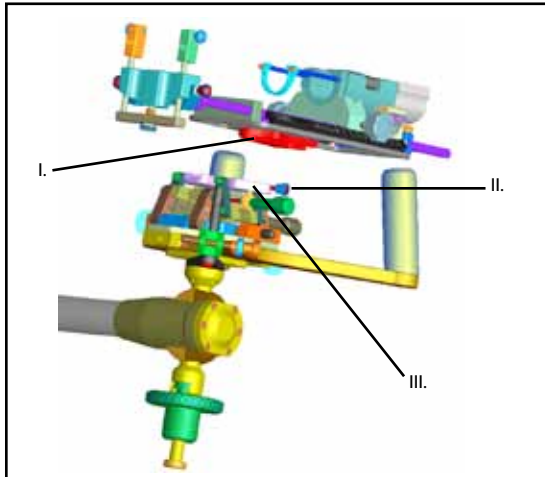
3G Ultrasound SimulView

- I. Grid platform
- II. Probe positioning locator
- III. Cradle latch
- IV. Cradle
- V. Grid rail locking knobs
- VI. Grid rail handle
- VII. Carriage scale
- VIII. Carriage scale marker
- IX. Carriage knob
- X. Carriage

SETTING UP STEPPER

- I. Slide stepper interface plate into stabilizer interface plate. Secure by tightening stepper quick connect knob.

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2. Rotate cradle out of carriage.
3. Open the cradle latch.
4. **For GE Healthcare ERB transducer:** Place and orient transducer into the cradle, sliding transducer forward until it stops.
For Aloka UST-672-5 / 7.5, 3G Ultrasound SimulView transducers: Place and orient transducer into the cradle.
For Acuson ER7B, Hitachi EUP-U533, Siemens Endo-P II transducers: Slide transducer back into the cradle until it stops.
For B-K Medical 8658 (8558) transducer: Slide cable end of transducer back into the cradle, ensuring locating feature fits into notch.
5. Close cradle latch and tighten knob.
6. Replace cradle with transducer back into stepper by rotating into carriage. Center transducer in rotation by feeling the center line detent during rotation.
7. Position carriage in a central location on carriage scale by turning carriage knob until carriage scale marker is at 60.
8. Loosen grid rail locking knobs and pull back on grid rail handle so grid platform is out of the way for initial transducer positioning.

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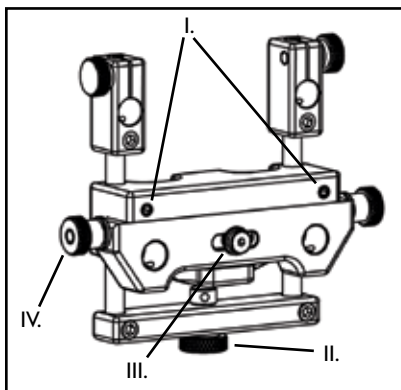
1. Grasp handles of FAM and loosen locking control knob on stabilizer.
For Micro-Touch: Adjust weight by turning weight adjustment knob on stabilizer to achieve desired feel to equipment.
2. Obtain satisfactory transducer position and starting image.
3. Install sterile drape and grid, position grid platform appropriately and secure with grid rail locking knobs.

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ADJUSTING GRID PLATFORM

1. Loosen Vertical Adjustment Screws with the provided hex wrench (3/32"). Align needle with the software by turning the Vertical Adjustment Knob until the vertical needle position is correctly aligned with grid software display. Tighten Vertical Adjustment Screws.
2. Loosen Horizontal Adjustment Screw with the provided hex wrench (3/32"). Align needle with the software by turning the Horizontal Adjustment Knob until the horizontal needle position is correctly aligned with grid software display. Tighten Horizontal Adjustment Screw.

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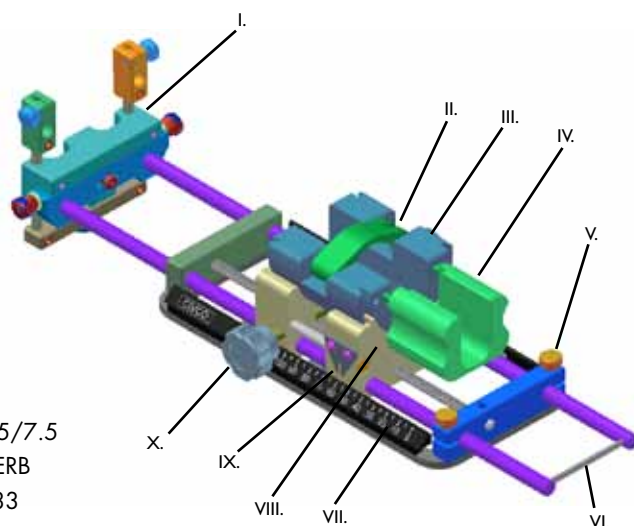
FUNCTIONAL CHECKS

- Prior to each use perform the following checks to ensure optimal stepper performance.
1. Rotate carriage knob to ensure carriage moves freely through full range of linear travel and performs stepping function.
 2. Rotate cradle to ensure movement through full range of rotational travel.

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UYARI

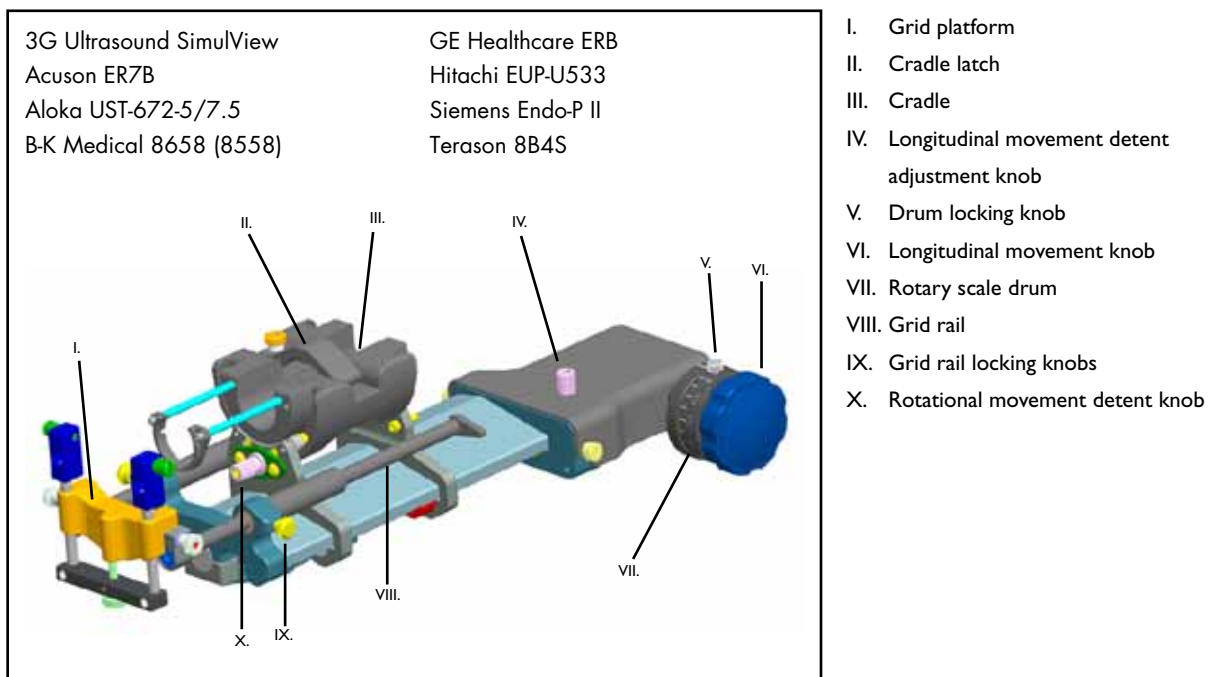
UYARI

Device Specification	Specification
Classification	Class I Type B
Operating Value	100-240 VAC
Current	250 mA -100 mA
Frequency	50 / 60 Hz
Protection Against Ingress of Water	Ordinary Protection
Power Supply	Model No: MW4012F
Transport Conditions	-29°C to 38°C
Storage Conditions	11°C to 29°C
Max. Relative Humidity	<85% RH

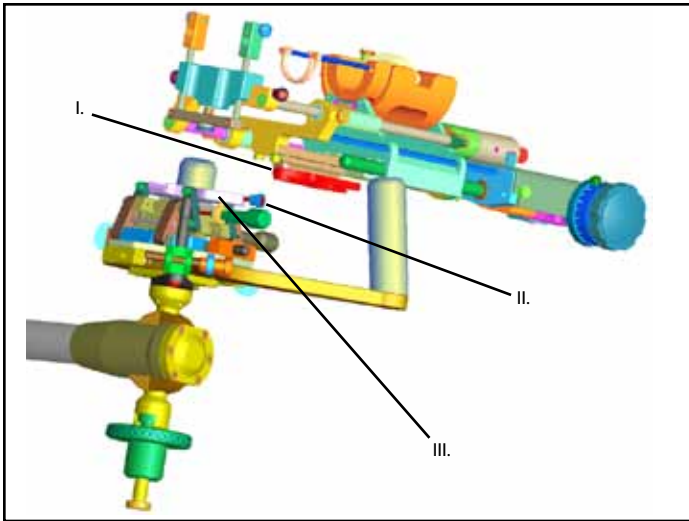
PRIOR TO USE

- The user is responsible for ensuring the application and use of the device does not compromise the patient contact rating of any equipment used in the vicinity of, or in conjunction with, the system.
- The use of accessory equipment and/or hardware not complying with the equivalent product safety and EMC requirements of this product may lead to a reduced level of safety and/or EMC performance of the resulting system. Consideration relating to the choice of accessory equipment used with this product shall include:
 - The use of the accessory in the patient vicinity.
 - Evidence the safety certification of the accessory has been performed in accordance with the appropriate IEC 60601-1 and/or IEC 60601-1-1 Harmonized National Standards.
 - Evidence the EMC certification of the accessory has been performed in accordance to the IEC 60601-1-2 Harmonized National Standards.
- Observe all safety precautions recommended by the accessory equipment manufacturer in the user documentation provided with the equipment.
- Keep this manual with device for reference as required.

SETTING UP STEPPER

I. Slide stepper interface plate into stabilizer interface plate. Secure by tightening stepper quick connect knob.

- I. Stepper interface plate
- II. Stepper quick connect knob
- III. Stabilizer interface plate

2. Open the cradle latch.

NOTE

• **For GE Healthcare ERB transducer:** Transducer handle must be removed.

3. Place and orient transducer into the cradle.

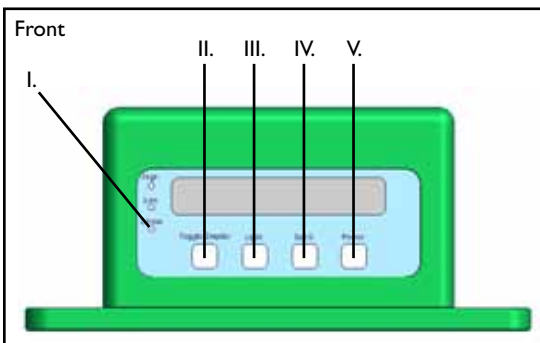
4. Close cradle latch and tighten knob.

5. Loosen grid rail locking knobs and pull back on grid rail so grid platform is out of the way for initial transducer positioning.

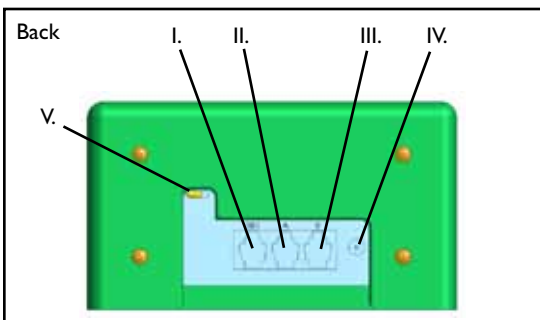
6. Center stepper using longitudinal movement knob.

NOTE

- Free longitudinal movement achieved by fully disengaging longitudinal movement detent knob by turning counter clockwise. Enable stepping function by tightening knob.
- Free rotation of cradle achieved by disabling rotational movement detent knob by turning clockwise. Enable rotational detent by tightening knob.

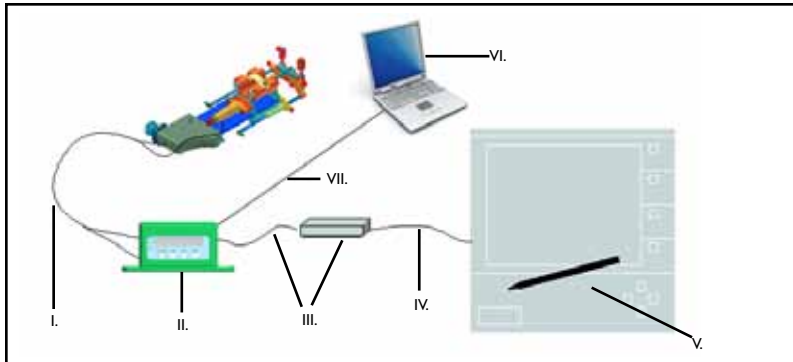
CONNECTING DIGITAL DISPLAY

- I. Linear Light - On = Linear Position / Off = Rotational Position.
- II. Toggle display button - Toggles display between rotational and linear display.
- III. Light button - Turns back light display on and off. Turn back light on in low light situations.
- IV. Set 0 button - Sets current display position to zero.
- V. Power button - Turns power on and off display.

- I. Serial interface port data output
- II. Encoder cable A connection port
- III. Encoder cable B connection port
- IV. Power supply port (12VDC)
- V. Programming toggle switch (no function for end user, CIVCO use only)

1. Connect encoder cable set to back of digital display. Connect power supply pack to the back of digital display. Connect the power cord between the power supply pack and the auxiliary power outlet on back of ultrasound system. Use power supply pack and power cord provided.
2. Optional: Connect serial cable to digital display and to UL 60950-I approved computer. Serial cable can be used to transmit information on position from digital display.

- I. Encoder cable set
- II. Digital Display
- III. Power supply pack (12VDC)
- IV. Power cord
- V. Ultrasound system auxiliary power outlet
- VI. Computer
- VII. Serial cable

ADJUSTING THE STEPPER AND STABILIZER

1. Grasp handles of FAM and loosen locking control knob on stabilizer.
- For Micro-Touch:** Adjust weight by turning weight adjustment knob on stabilizer to achieve desired feel to equipment.
2. Obtain satisfactory transducer position and starting image.
3. Install sterile drape and grid, position grid platform appropriately and secure with grid rail locking knobs.
4. To manually set longitudinal scale, loosen drum locking knob and position rotary scale drum at zero. Tighten drum locking knob to secure.
5. To digitally set rotational scale, place cradle in center position and reset rotational position display to zero, refer to front digital display illustration.
6. Press the toggle display button to change display to longitudinal position.
7. To digitally set longitudinal scale, set digital display at zero by pressing the Set 0 button.

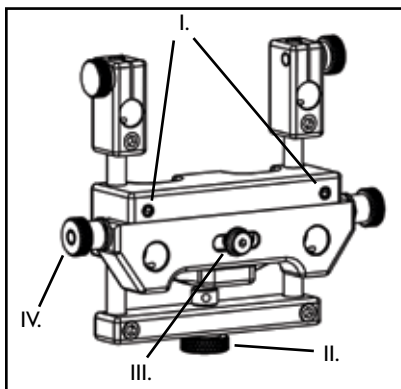
WARNING

- Before using assembled equipment, verify needle position is correctly aligned with grid software display. If needle position is not correctly aligned with grid software display, adjust grid platform.

ADJUSTING GRID PLATFORM

1. Loosen vertical adjustment screws with the provided hex wrench (3/32"). Align needle with the software by turning the vertical adjustment knob until the vertical needle position is correctly aligned with grid software display. Tighten vertical adjustment screws.
2. Loosen horizontal adjustment screw with the provided hex wrench (3/32"). Align needle with the software by turning the horizontal adjustment knob until the horizontal needle position is correctly aligned with grid software display. Tighten horizontal adjustment screw.

- I. Vertical Adjustment Screws
- II. Vertical Adjustment Knob
- III. Horizontal Adjustment Screw
- IV. Horizontal Adjustment Knobs

DEVICE INSPECTION

1. Stepper should hold transducer securely, allow smooth rotation and precise movement along longitudinal axis of the transducer.

FUNCTIONAL CHECKS

- Prior to each use perform the following checks to ensure optimal stepper performance.

1. Rotate longitudinal movement knob to ensure carriage moves freely through full range of linear travel and performs stepping function.
2. Rotate cradle to ensure movement through full range of rotational travel.

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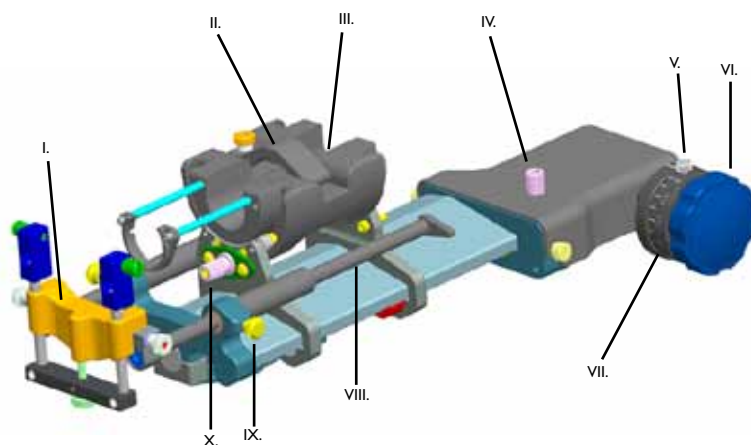
3G Ultrasound SimulView
Acuson ER7B
Aloka UST-672-5/7.5
B-K Medical 8658 (8558)

GE Healthcare ERB
Hitachi EUP-U533
Siemens Endo-P II
Terason 8B4S

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WARNING**REPROCESSING**

- Users of this product have an obligation and responsibility to provide the highest degree of infection control to patients, co-workers and themselves. To avoid cross-contamination, follow infection control policies established by your facility,
- Do not soak or immerse the EXII model.
- Do not place in a mechanical washer.
- Do not use ethylene oxide or autoclave to sterilize.

1. Manually soak non-electrical parts of stepper with mild detergent, such as ENZOL® Enzymatic Detergent (Johnson & Johnson) following the manufacturers recommendation for use.
2. Disinfect all surfaces with common germicidal or antiseptic wipe, such as alcohol or hydrogen peroxide.
3. Wipe dry with lint-free cloth or allow to air dry prior to reassembly and use.

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WARNING**MAINTENANCE**

- Device should be stored to keep all components together and secure.
- If device does not function per instructions, contact CIVCO.

NOTE**SHIPPING AND TRANSPORTATION**

- All product to be returned must be in its original packaging. Contact CIVCO for further instructions as needed.

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NOT

Drape

INTENDED USE | • The drape is to be used as an instrument cover for medical devices and helps to prevent transfer of microorganisms, body fluids and particulate material.

WARNING

- In the USA, federal law restricts device to sale or use by or on order of a physician.
- Before use, you should be trained in ultrasonography. For instructions on the use of your transducer, see your system's manual.
- Do not use if integrity of packaging is violated or if expiration date has passed.

APPLYING DRAPE

NOTE | • To ensure sterile procedure, drape unit after transducer is accurately positioned in rectum and radiation planning is completed.

Taped Drape

1. Using sterile technique, remove drape from package.
2. Unfold drape with edge markings up. Tear tab on cut side.
3. Place drape over stepper with cut side open toward patient.
4. Secure with tape strips under transducer.
5. Align paired holes in drape with holes in grid platform.
6. Securing grid to grid platform will secure drape in position.

Gathered Drape

1. Using sterile technique, remove drape from package.
2. Unfold drape and locate circular stickers on top corner.
3. Stretch bottom banded portion of drape and place over stepper with sticker side toward patient.
4. Align cross-hair punches within stickers on drape with holes in grid platform.
5. Secure grid to grid platform through cross-hair punches in drape.

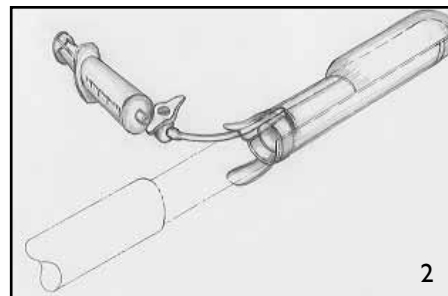
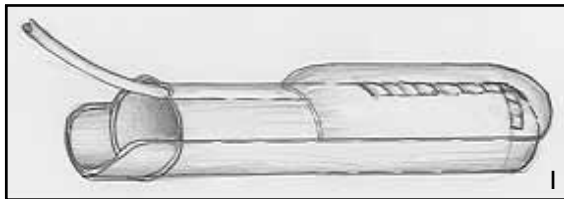
DISPOSAL

WARNING

- Dispose of single-use components as infectious waste.

Endocavity Balloon

INTENDED USE | • The transducer standoff device can be used to perform scanning and needle guided procedures for body surface and endocavity diagnostic ultrasound.

WARNING

- In the USA, federal law restricts device to sale or use by or on order of a physician.
- Before use, you should be trained in ultrasonography. For instructions on the use of your transducer, see your system's manual.

Applying Balloon to Transducer

1. Place an appropriate amount of gel inside endocavity balloon and/or on transducer face to aid in scan quality, application and removal of balloon.
2. Insert transducer into balloon with fill tube oriented in top-center of transducer. Pull balloon tightly over transducer face to remove wrinkles and air bubbles.
3. Visually orient balloon so lateral seams are symmetrical distance from longitudinal array.
4. Tape balloon securely making sure to attach around plastic tabs (Tape not included).
5. Fill syringe with 30cc of saline solution. Fill balloon through luer lock port holding probe tip down; fill and aspirate with 10-20cc until all air is removed. Do not over inflate balloon.
6. Deflate balloon leaving 25-30cc syringe connected.

Disposal

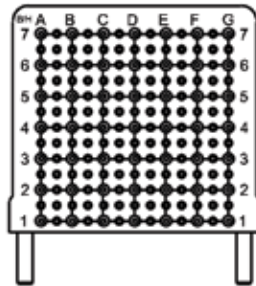
WARNING

- Dispose of single-use components as infectious waste.

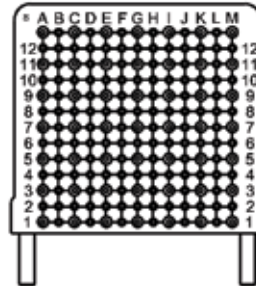
Grid

B-K Medical 8551, 8658 (8558) transducers
Hitachi EUP-U533 transducer

17/18GA

Acuson ER7B transducer
Aloka UST-672-5 / 7.5 transducer
GE Healthcare ERB transducer
Siemens Endo-P II transducer

17/18GA

INTENDED USE

- Provides physicians with a series of precisely placed needle guide channels used to facilitate accurate placement and insertion in relation to the ultrasound imaging probe.

WARNING

- In the USA, federal law restricts device to sale or use by or on order of a physician.
- Before use, you should be trained in brachytherapy. For instructions on preparation, safe operation and care of your transducer and brachytherapy system, see your system's manual.
- Do not use if integrity of packaging is violated or if expiration date has passed.
- Do not use grid if it appears to be cracked or broken.
- Grid cannot be resterilized after use.
- Ensure graphics on ultrasound monitor match graphics on grid.
- Take care when inserting needles into grid to avoid breaking needle.

ASSEMBLING GRID

- Using sterile technique, remove grid from package.
- Select appropriate graphics on grid to match ultrasound monitor.
- Tighten locking knobs.

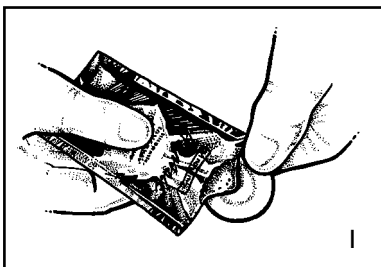
DISPOSAL

- Dispose of single-use components as infectious waste.

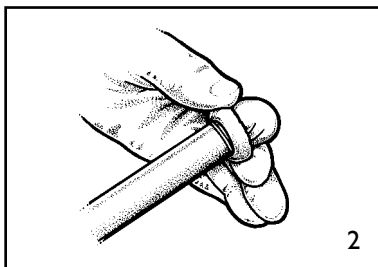
WARNING

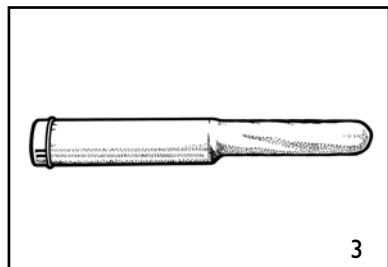
NeoGuard™ Latex-Free Transducer Cover

1

2

3

INTENDED USE

- The cover allows use of the transducer in scanning and needle guided procedures for body surface, endocavity and intra-operative diagnostic ultrasound, while helping to prevent transfer of microorganisms, body fluids and particulate material to the patient and health care worker during reuse of the transducer.

WARNING

- In the USA, federal law restricts device to sale or use by or on order of a physician.
- Before use, you should be trained in ultrasonography. For instructions on the use of your transducer, see your system's manual.
- Dispose of single-use components as infectious waste.
- See your system's manual for reprocessing transducer between use.

WARNING

Covering the Transducer and Bracket

- Use only water-soluble agents or gels. Petroleum or mineral oil based materials may harm cover.
- Disposable components are packaged sterile and are single-use only. Do not use if integrity of packaging is violated or if expiration date has passed.
- For illustration purposes only, transducer may be shown without a transducer cover. Always place a cover over transducer to protect patients and users from cross-contamination.

- Place an appropriate amount of gel inside cover and/or on transducer face. Poor imaging may result if no gel is used.
- Insert transducer into cover making sure to use proper sterile technique. Pull cover tightly over transducer face to remove wrinkles and air bubbles, taking care to avoid puncturing cover.
- Secure cover with enclosed bands (if applicable).
- Inspect cover to ensure there are no holes or tears.

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Corporate Headquarters: 102 First Street South, Kalona, IA 52247-9589 USA
319.656.4447, 800.445.6741 | fax: 319.656.4451, 877.329.2482
www.civco.com

COPYRIGHT © 2011 ALL RIGHTS RESERVED. CIVCO IS A REGISTERED TRADEMARK OF CIVCO MEDICAL SOLUTIONS. ACCUCARE, MICRO-TOUCH AND NEOGARD ARE TRADEMARKS OF CIVCO MEDICAL SOLUTIONS. ALL OTHER TRADEMARKS ARE PROPERTY OF THEIR RESPECTIVE OWNERS. ALL PRODUCTS MAY NOT BE LICENSED IN ACCORDANCE WITH CANADIAN LAW. PRINTED IN USA. US PATENTS 5931786, 6036632, 6248101.

043-658H

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MPS GmbH